

READ ME FIRST

# QUICK START GUIDE

for illumia®Pro3


Labsphere

This quick-start guide provides the fundamental steps required to set up, connect, and power up the Labsphere illumiaPro3 System.

Match the black hardware component identification numbers with their corresponding **red** connection diagram numbers on page 2.

For complete information, refer to the full Labsphere illumiaPro3 hardware manual AQ-41000-050 and Integral Software manual AQ-81000-100 included on the Labsphere USB flashdrive. (item 8)

**ITEMS MAY BE IN SEPARATE BOXES WITHIN SHIPPING CRATE.  
ENSURE ALL ITEMS ARE PRESENT.**

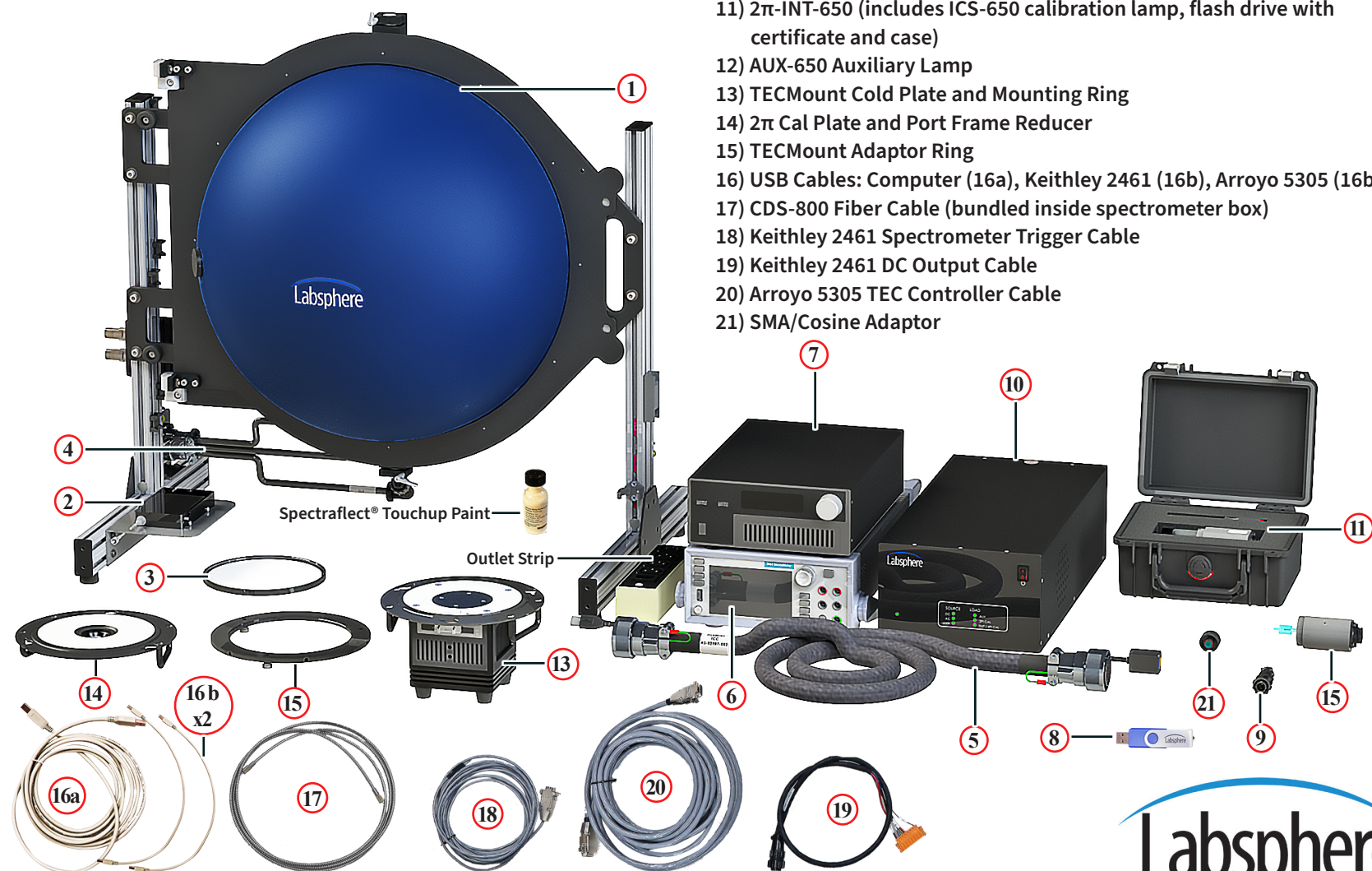


## HARDWARE COMPONENTS

- 1) Integrating Sphere Assembly
- 2) CDS-800 Spectrometer
- 3) 2 $\pi$  Port Frame Reducer (for optional Forward Flux Standard)
- 4) Lamp Cable Bundle and Junction Connector (pre-installed on #1 Sphere)
- 5) Interconnect Cable
- 6) Keithley 2461 SourceMeter®
- 7) Arroyo 5305 TEC Controller
- 8) Flash Drive: System documentation and Integral® software
- 9) ICM-500 Shorting Plug (install into #10 ICM-500)
- 10) Labsphere ICM-500
- 11) 2 $\pi$ -INT-650 (includes ICS-650 calibration lamp, flash drive with certificate and case)
- 12) AUX-650 Auxiliary Lamp
- 13) TECMount Cold Plate and Mounting Ring
- 14) 2 $\pi$  Cal Plate and Port Frame Reducer
- 15) TECMount Adaptor Ring
- 16) USB Cables: Computer (16a), Keithley 2461 (16b), Arroyo 5305 (16b)
- 17) CDS-800 Fiber Cable (bundled inside spectrometer box)
- 18) Keithley 2461 Spectrometer Trigger Cable
- 19) Keithley 2461 DC Output Cable
- 20) Arroyo 5305 TEC Controller Cable
- 21) SMA/Cosine Adaptor

Labsphere

SYSTEM SETUP

- 1. Remove all components from the shipping crate and identify them with image callouts on Page 1 and number stickers on the individual package items.
- 2. Remove the shipping lockbolt above the clasp handle.
- 3. Place the electronic components (items 6, 7 and 10) to be within 14ft (4.3M) or less of the junction box on the sphere.
- 4. Put on nitrile or latex gloves, then connect the system's components and cables using the item numbers in the Hardware Components list on the front page as reference in the Connection Diagram below.  
**Note:** Refer to the illumiaPro3 hardware manual AQ-41000-050 on the flash drive [item 8] for detailed assembly information.
- 5. Turn on the host computer and wait for Windows to load, then turn on the power switch on the front panel of each electronic component.
- 6. Insert the supplied Labsphere USB flash drive into a USB socket in the host computer and run the Integral Installer.  
**Note:** Refer to the Integral software manual AQ-81000-100 on the flash drive (item 8) for detailed information.

 Components may be contained in multiple boxes within the crate.

