

READ ME FIRST

QUICK START GUIDE

Integral Software
Installation Quick-Start




Labsphere Integral® Software Installation Quick-Start

This document covers the installation of the Labsphere Integral software. For complete information, refer to the Labsphere illumia hardware manual and the Integral software manual AQ-81000-100 included in the Labsphere USB flash drive.

SOFTWARE INSTALLATION

1. Connect and power on the system components as defined in the illumia hardware manual.
2. Insert the Labsphere flash drive into a USB port in the host Windows 10 or higher computer.
3. Navigate down to the compressed “IntegralInstall” folder and right-click on it to extract the installer files to the computer desktop.
4. Right click on the “IntegralInstall64Bit” Windows batch file and select “Run as administrator”.
5. The Python 2.7 (beginning) and LSAPI installer (end) run separately. Popup windows will close automatically when they are complete. Choose (D) when prompted for the file (F) or directory (D).
6. Install the Chrome internet browser if it is not on the host computer.
32-bit Chrome is included in the <\Integral Install\Assets\> folder (ChromeSetup32bit.exe).
7. The Integral (URL) shortcut will appear on the desktop. This shortcut will open in the default browser (Integral requires Chrome). Update the Windows default browser if necessary, in [System Settings -> Default apps].

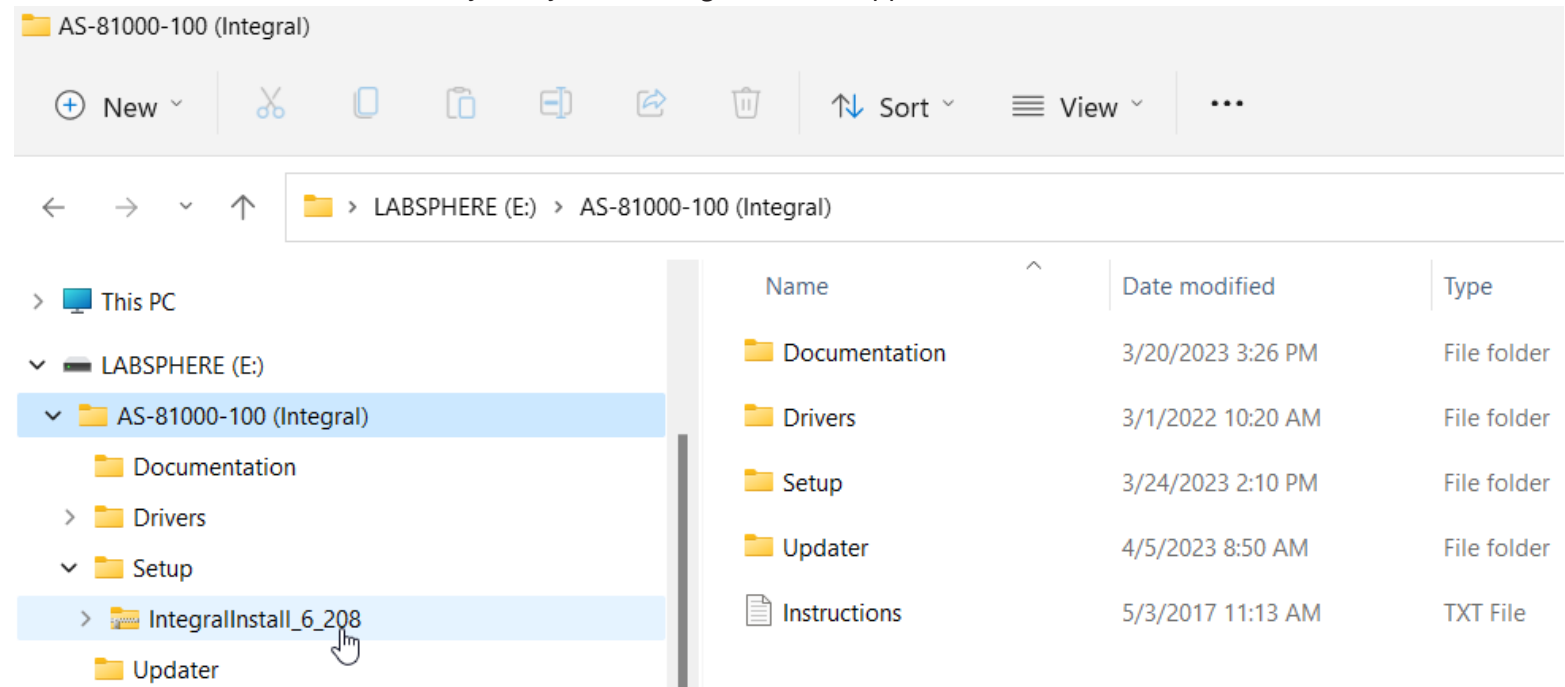

OPEN INTEGRAL®

1. Open the Chrome browser when Windows has loaded. Integral can also be manually launched via the Chrome URL=(localhost:5961). For local control, type “localhost:5961 in the address line to access the Integral software. For remote control, type the IP address, colon, then the port on any other computer that can access the network. For example, 192.168.1.10:5961.
2. Log in by using “**admin**” as the name and password.





www.labsphere.com | +1 603.927.4266

© 2023 Labsphere, Inc. Trademark Names® are Registered Trademarks of Labsphere, Inc. All Rights Reserved